




SONICA® 3200L

6 litres

Ultrasonic cleaning tank available in version M, MH, ETH and EP. Suitable for use in chemical laboratories when cleaning particularly long objects such as pipettes, cuvettes, glassware, capillary tubes, electrodes and endoscopes for medical use.

Tank size

500x140x100mm

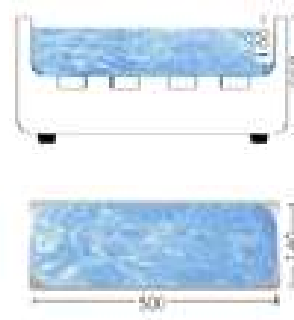

General informations

- Average mains power when not using heater 180W
- Average mains power when using heater 355W
- HF Power 720W p-p.
- Ultrasound frequency 45Khz \pm 5%

- Standard drain
- Number of transducers 4
- Capacity 6lts.
- Case and tank made in stainless steel AISI 304

Control panel



VERSION M
The simplest version in the SONICA® range is equipped with a timer which can be set for to 15 minutes. It can also be set to manual operation.



VERSION ETH
Cleaning times and temperatures controlled electronically by a microprocessor.
Waterproof control panel.



VERSION MH
This is equipped with a timer which can be set for up to 15 minutes and can also be set to manual operation. The thermostat for the heater is set at 60°C.



VERSION EP
The most technologically advanced of the SONICA® range.

- New Sweep system generator.
- Modern digital display control panel windows.
- Degas time setting.
- Cleaning times from 0 to 99 minutes.
- Tank liquid temperature setting up to 70°C.
- Memory for 9 cleaning programmes.
- Jet program function to start quickly and easily cleaning cycle.

Accessories

- Stainless steel holder basket
- Small stainless steel immersion basket
- Glass beaker
- Plastic beaker
- Stainless steel lid

Liquid detergents

SONICA® UG Concentrated detergent for general use. Bottle Kg1.

SONICA® SF Concentrated detergent for removing heavy grease from very dirty spots. Bottle Kg1.

SONICA® AF Concentrated detergent for removing limestone deposit stains and large quantities of oxides and rust. Bottle Kg1.

SONICA® AC Concentrated detergent for removing limestone traces and small amounts of oxides and rust. Bottle Kg1.

SONICA® RC Ready-to-use detergent for removing cement deposits from dental instruments. Bottle Kg1.

SONICA® RA Ready-to-use detergent for removing alginate deposits from impression trays in dental application. Bottle Kg1.

SONICA® RG Special neutral detergent for removing plaster deposits. Bottle Kg1.

